## Ultra-Compact Digital Mixer

### English

Instruction Manual  
pages 3-8

### Français

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### Espanol

Manual de instrucciones  
pages 15-20

### Deutsch

Bedienungsanleitung  
pages 21-26

# Contents

|                                              |     |
|----------------------------------------------|-----|
| GENERAL INSTRUCTIONS .....                   | 4   |
| SAFETY INSTRUCTIONS .....                    | 4-5 |
| APPROVALS .....                              | 5   |
| TECHNICAL DATA.....                          | 5   |
| MIXER DESCRIPTION.....                       | 6   |
| INSTALLATION .....                           | 6   |
| ENVIRONMENTAL OPERATING CONDITIONS .....     | 7   |
| GENERAL OPERATING INSTRUCTIONS .....         | 7   |
| CLEANING, SERVICE, TRANSPORT & STORAGE ..... | 7-8 |
| WARRANTY AND PRODUCT LIABILITY .....         | 8   |

## PACKING LIST

- |                                  |                            |
|----------------------------------|----------------------------|
| • Mixer                          | • Impeller shaft           |
| • Power cord                     | • (2) Impeller blades      |
| • Stand base                     | • (2) Impeller lock screws |
| • Stand support rod and lock nut | • Instruction manual       |

## GENERAL INSTRUCTIONS

Read all safety instructions before operating.

Unpack the mixer carefully. Ensure the mixer is not damaged and verify the contents against the packing list. If the mixer is damaged or the contents do not correspond with the packing list, notify your supplier immediately.

Read the manual carefully! Make sure that all users read the manual thoroughly before operating the mixer. For further questions or inquiries, please contact:

Cole-Parmer

625 East Bunker Ct.

Vernon Hills, IL 60061-1844

Tel: 800-323-4340 (within the U.S. and Canada),

847-549-7600 (outside the U.S. and Canada)

Web site: [www.coleparmer.com](http://www.coleparmer.com)

E-mail: [info@coleparmer.com](mailto:info@coleparmer.com)

## SAFETY INSTRUCTIONS

Each operator should read all instruction before using this device.

The mixer must be securely fixed to a stable support such as the provided stand. If other stands are used, the unit must be secured such that it will not fall if tipped.

**CAUTION:** This is not an explosion-proof mixer. Do not use in a highly flammable or explosive environment.

Spinning paddles or impellers can cause severe personal injuries. Operators must use extreme care and good judgment when mixing at any speed.

All mixing paddles and impellers must be in good condition with straight shafts. If the mixer vibrates at high speeds, check the paddle shaft for damage, and repair or replace if needed.

Extreme care must be taken when mixing chemicals to ensure that no chemicals are splashed outside the mixing vessel. Care must be taken when changing to faster mixing speeds.

Ensure that the mixing impeller does not contact the containment vessel.

Do not operate while standing in water. Keep the unit dry and do not immerse any part, except the mixing paddle into any liquids. Protect unit and operator from splashing.

Ensure that no loose clothing, jewelry, or hair can become entangled in any rotating parts. A fast spinning chuck can cause injury to operator.

Power can be interrupted to the mixer by pressing the power button. In an emergency, disconnect the main power cord.

## APPROVALS

### North American

Constructed in accordance with CAN/CSA-C22.2 (1010),  
FCC 15/ICES-003/EN 55011

Tested to the same standards as IEC61010.

### European

Constructed in accordance with EMC/EMI standard EN61326-1:2006 (EN55011, conducted and radiated emissions, EN61000-3-2, EN61000-4-2, EN61000-4-3, EN61000-4-4, EN61000-4-5, EN61000-4-6)


RoHS and WEEE compliant.

## TECHNICAL DATA

|                                              |                                                              |
|----------------------------------------------|--------------------------------------------------------------|
| <b>Model Number</b>                          | <b>50006-00</b>                                              |
| Speed range                                  | 50 to 2500 rpm                                               |
| Maximum torque                               | 10 Ncm                                                       |
| Mechanical output power                      | 1/100 hp, 7 watts                                            |
| Input to power supply                        | AC 100 to 240 volts, 50/60 Hz (must use supplied power cord) |
| Electrical output from power supply to mixer | DC 24 V, 0.5 A                                               |
| Weight - Mixer only                          | 0.9 lbs (.4 kgs)                                             |
| Weight - Total package                       | 3.1 lbs (1.42 kgs)                                           |
| Maximum volume                               | 2 liters, water-like viscosity                               |
| Molded housing                               | Chemical-resistant, high-impact nylon. Recyclable (PA).      |
| Keypad                                       | Waterproof, chemical-resistant polyester                     |
| Power cord length                            | 6 ft                                                         |
| Keyless chuck                                | Fits 6 mm shaft                                              |
| Led display accuracy                         | ±10 rpm, 4 digit                                             |

Technical specifications subject to change without notice.

## MIXER DESCRIPTION

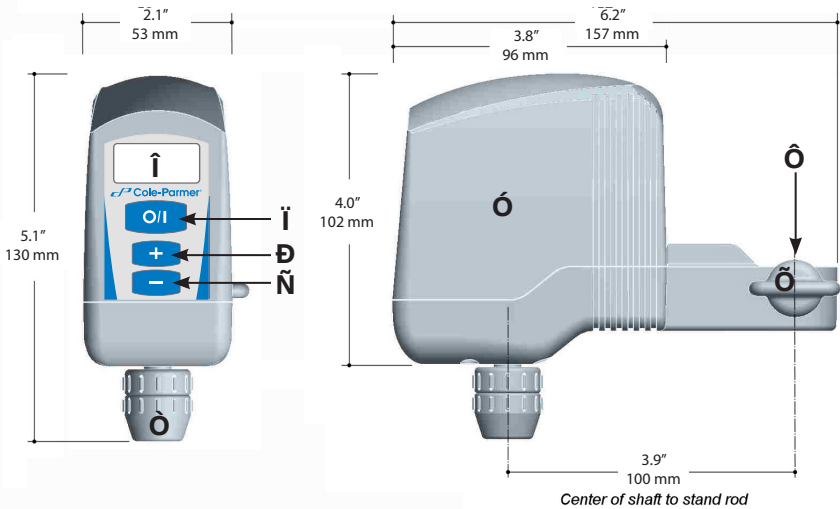- Î Four-digit LED display
- İ Button - power on/off
- Đ Button - increase rpm
- Ñ Button - decrease rpm

- Ò Keyless chuck
- Ó Housing
- Ô Clamp opening
- Õ Clamp knob

## INSTALLATION

### Stand Assembly

Remove lock nut. Insert the threaded end of the stand rod into the stand base. Turn rod clockwise to tighten. Secure with lock nut.

### Mixer Installation

Slide the clamp opening over the top of the stand rod. Lower to desired height. Turn clamp knob to securely tighten the mixer onto the stand rod. Connect supplied power cord to mixer. Connect power cord to a grounded electrical outlet.

### Impeller Assembly

Add impeller blade(s) to impeller shaft. Tighten set screw to secure.

### Impeller Installation

Turn lower half of chuck to open. Do not remove completely. Insert impeller end into chuck opening until it stops. Turn lower half of chuck (finger tight) to close chuck and secure impeller.

## ENVIRONMENTAL OPERATING CONDITIONS

The mixer must operate in the following conditions:

- Indoors
- Altitudes not exceeding 2000 meters
- Temperature from 41 to 104°F (5 to 40°C)
- Maximum relative humidity 80% for temperatures up to 87°F (31°C) decreasing linearly to 50% relative humidity at 104°F (40°C)
- Pollution degree 2 in accordance with IEC 664 "Normally only non-conductive pollution occurs. Occasionally, however, temporary conductivity caused by condensation may be expected."

## GENERAL OPERATING INSTRUCTIONS

The mixer must be securely mounted as specified in the installation instructions. Ensure the mixing shaft and blade are securely attached. The mixer is now ready for use to mix liquids.

Button Operation Functions:

Power Button - Turns the mixer motor on and off.  
When initially turned on, mixer will spin at 50 rpm.

Increase Button - Increases the speed.

Decrease Button - Decreases the speed.

If the load is too great, an overload condition may occur and the mixer will automatically power off. To resume operation, reduce the mixer load and press the power button.

## CLEANING, SERVICE, TRANSPORT & STORAGE

### Cleaning

Using a mild soap and water solution, wipe the housing and control panel of the mixer with a damp cloth. Do not use chlorine bleach, chlorine-based cleanser, abrasives, ammonia, steel wool, or scouring pads with metal content or similar harsh solvents or abrasives. These may damage the surface of the instrument.

### Service

Do not attempt to service or repair mixer. **If the mixer housing is opened, the warranty becomes void.** Contact Supplier for return authorization and return instructions.

### Transport and Storage

Clean mixer so that it is free from any materials which may be harmful. Provide an MSDS or other safety data sheets for all chemicals used in conjunction with this mixer where appropriate. Place the stirrer and its parts into the original packing to protect and prevent damage during transport. Seal the original packing with packing tape. Store the packed stirrer in a dry place.

**CAUTION:** Failure to clean, service, and handle the mixer as outlined above can cause it to fail which could impair the safety of the user.

## **WARRANTY AND PRODUCT LIABILITY**

### Warranty

Cole-Parmer warrants the herein described product of its manufacture for three years from receipt of the equipment by the purchaser, against defects in material and workmanship. For warranty registration fill out the supplied warranty registration card.

This limited warranty covers parts and labor. **If the mixer housing is opened, the warranty becomes void.**

This warranty may only be altered by a specifically published amendment. No individual has authorization to alter the provisions of this warranty policy or its amendments. Expendable items are not covered by this warranty. Damages due to corrosion or accidental breakage are not covered.

### Liability

Under no circumstances shall Cole-Parmer be liable for indirect, consequential, or special damages of any kind.

Cole-Parmer reserves the right to make technical changes without prior notice.

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